



Scooterland
MOBILITY
Gold Coast

07 5539 2141

VECTOR REHAB CHAIR

Multi-Function Rehab Power Wheelchair

- Center of gravity changes according to sliding seat. (For safety reasons, centre gravity of end-user will always stay above the middle wheels)
- Mechanical shear reduction device for injury prevention
- Armrest adjusted automatically when reclined for comfortable seating
- Adjustable seat width and seat depth (16"/18"/20") for customized fitting of individual user
- Dual in-line motors for enhanced efficiency, torque, range and performance
- Optional seat elevator and power footplate or elevating leg rest

Optional PG R-NET multi-function controller



165 Currumburra Road Ashmore QLD 4214
Unit 6 (Showroom), Unit 2 (Workshop)
sales@scooterland.com.au

www.scooterland.com.au

NDIS




Scooterland
MOBILITY
 Gold Coast

07 5539 2141

**VECTOR
 REHAB CHAIR**



Specification - Vector Rehab Chair	
Combination Options;	Combination Descriptions;
MRT	Manual Recline + Power Tilt
PRT	Power Recline + Power Tilt
SET	Seat Elevator + Power Tilt
ERT	PRT (Power Recline + Power Tilt) + Seat Elevator

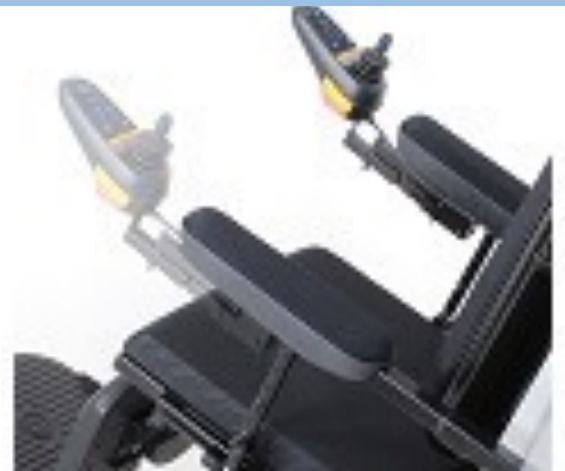

Specification - Vector Rehab Chair	
Weight Capacity	182kg / 400lb
Max. Speed	8kph / 5mph
Range	38km / 24mi.
Turn Radius	68cm / 27"
Ground Clearance	8cm / 3"
Motor	DC 24V, 200W (Rating) / 575W (Max)
Controller	PG VR2 60A / R-Net 60A
Battery	22NF / 12V 60A/Hr x 2pcs
Charger	5A off-board
Gradient	7.5°
Caster	Front: 8" & Rear: 8" (200mm x 50mm)
Drive Wheel	14"x 3" (3.00-8)
Brake	Intelligent, Regenerative, Electromagnetic Brakes

